SUN NUCLEAR
A MIRION MEDICAL COMPANY

Multi-Energy CT Phantom

Comprehensive Testing, Tissue Equivalence

- Evaluate spectral CT performance with confidence: quantitative QA for dual-energy and photon-counting
- Features 28 inserts representing different dimensions and concentrations of iodine, calcium, blood, adipose and other materials of particular interest to Multi-Energy CT (MECT)
- Expanded range of HE Iodine and HE Calcium




Multi-Energy CT scanners have enabled improved clinical differentiations, such as distinguishing blood from calcification and calcification from iodinated contrast.[1,2] They can also create virtual mono-energetic images for clinical evaluation. However, the ability to achieve these benefits depends not only on one's equipment, but on the protocols used.

Sun Nuclear's Multi-Energy CT Phantom enables robust evaluation of scanner performance.

- Test material discrimination using solid rods representing iodine, calcium, blood, adipose, and more
- Ensure the efficacy of clinical protocols for multi-energy analysis
- Verify the quantitative accuracy of multi-energy scans
- Compare the consistency and stability across different scanners
- Check for artifacts in an extended field-of-view
- Test in both head (20 cm) and body (40cm x 30 cm) configurations
- Enable automated analysis with patented rod marker technology

Ensure Accuracy

Scanner accuracy can vary based upon scanner hardware, the dual energy post-processing, and the mAs used. Without an appropriate phantom, neither scanner accuracy nor variability are well known. Use of our calibrated high-Z inserts enables such quantification. Additionally, protocols that appropriately balance patient dose with system performance can be identified.

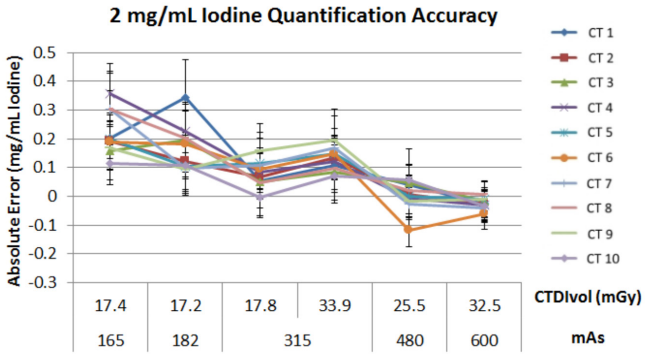

2mg/mL iodine rod quantification accuracy by mAs for 10 fast-kVp switching Multi-Energy CT scanners. Accuracy averaged over a 1 year period.

sunnuclear.com

©2026 Mirion Technologies, Inc. or its affiliates. All rights reserved. Sun Nuclear, the Sun Nuclear logo, and other trade names of Sun Nuclear products listed herein are registered trademarks or trademarks of Mirion Technologies, Inc. or its affiliates in the United States and other countries. Third party trademarks mentioned are the property of their respective owners. All Rights Reserved. All data used is best available at time of publication. Data is subject to change without notice.


MM-0101 - 11/2026

Automated analysis


The Multi-Energy CT Phantom supports automated analysis using patented rod marker technology. Each insert is tagged with a pattern, making it uniquely identifiable in a CT scan. Upcoming software will leverage these identifiers to evaluate results quickly and with minimal user interaction.

Enhance confidence in your virtual mono-energetic images

Monochromatic HU numbers have been shown to vary between scanners.[3] Moreover, the performance of multi-energy algorithms can be compromised by insufficient mAs. By utilizing Tissue Mimicking Material that replicates expected HU dependencies from 40-200 keV, the Multi-Energy CT Phantom lets you quantify these effects and define effective operating parameters.

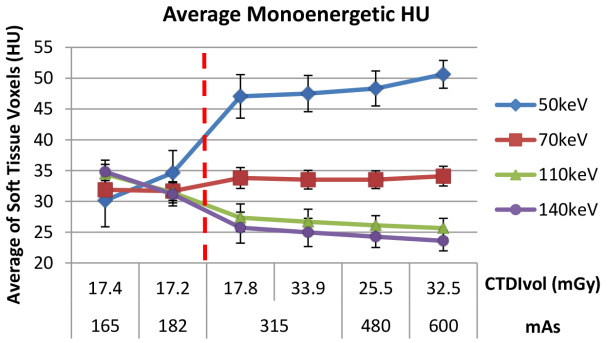

Average soft tissue monoenergetic HU for 50, 70, 110, and 140 keV reconstructions versus mAs, averaged over a 1 year period. Protocols left of the red dashed line were insufficient to provide reliable HU values.

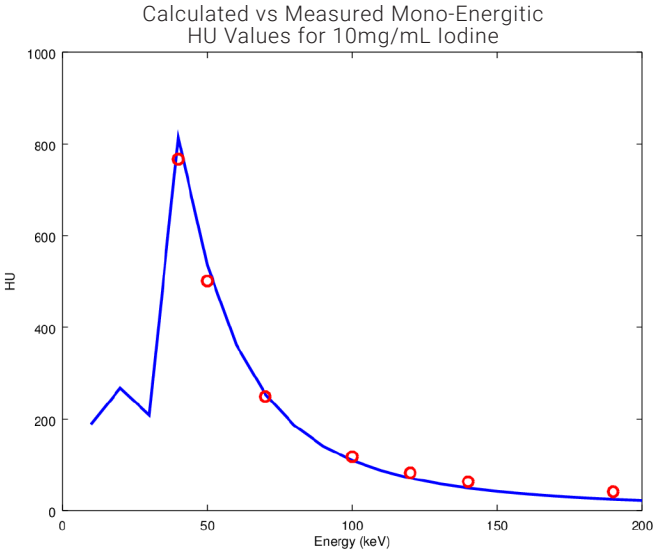

HU values of iodinated rod for mono-energetic reconstructions. Calculated values based on material compositions and NIST values (blue curve) vs HU values from mono-energetic reconstructions (red circles).

Evaluate an extended field-of-view

The ACR Quality Control Manual recommends checking for artifacts in a larger phantom on a weekly or monthly basis. The 40 cm extended field size of the Multi-Energy CT Phantom enables this artifact check to be performed concurrently with other evaluations, efficiently fitting into your workflow.

Specifications

In-plane Dimensions:	40.0 cm (15.7 in) x 30.0 cm (11.8 in)
Depth:	16.5 cm (6.3 in)
Diameter of Removable Head Section:	20.0 cm (7.87 in)
Material:	HE CT Solid Water®
Interchangeable Inserts:	27 solid inserts plus 1 true water container, each tagged with a CT-visible rod identification code
8 HE Iodine Inserts with Variable Concentrations:	Concentrations of 0.2, 0.5, 1.0, 2.0, 5.0, 10.0, 15.0, and 20.0 mg/mL
3 Iodine Inserts with Variable Diameters:	5.0 mg/mL concentration at diameters of 2.0, 5.0, and 10.0 mm
8 HE Calcium Inserts with Variable Concentrations:	Concentrations of 0, 5, 10, 20, 50, 100, 200, and 300 mg/mL
3 Blood [iron] Inserts:	Blood-mimicking material at relative electron densities of 1.03, 1.07, and 1.10
2 Blood [iron] with Iodine Inserts:	Blood-mimicking material plus iodine at 2.0 and 4.0 mg/mL
3 Tissue-Mimicking Inserts:	High-Equivalency Brain, High-Equivalency Adipose, High-Equivalency CT Solid Water
Weight:	15.5 kg (34.1 lbs)
Wheeled Case:	Included
Stand:	Included

DATA PROVIDED BY UT MD ANDERSON

[1] Nute JL, Jacobsen MC, Chandler A, Cody DD, Schellinghout D, Dual-Energy Computed Tomography for the Characterization of Intracranial Hemorrhage and Calcification: A Systematic Approach in a Phantom System. Invest Radiol. 2016; Jul 1

[2] Knoss N, Hoffman B, Krauss B, et al. Dual energy computed tomography of lung nodules: Differentiation of iodine and calcium in artificial pulmonary nodules in vitro. Eur J Radiology. 2011; 80(3): E516-519

[3] Mileto A, Barina A, Marin D, Stinnett S, Choudhury K, Wilson J, Nelson R Virtual monochromatic images from dual-energy multidetector CT: Variance in CT numbers from the same lesion between single-source projection-based and dual-source image-based implementations Radiology 2016 (in press)